


## Settings




Language

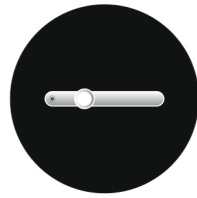

Brightness

## Instruction of remote controller

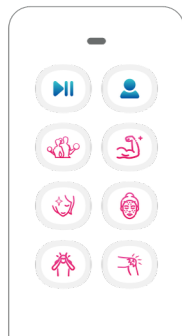

 Pause/Play

 Customized mode

 Fitness mode

 Recovery mode

 Skin care mode

 Acne mode

 Pain relief mode

 Inflammation mode

Before using the remote controller, please check if battery is included.



Battery normally not included for shipping reason

## Warranty Clauses

- 1) Super series are warranted against defects in materials and workmanship for a period of ONE(1) YEAR from the purchase date with free replacement.
- 2) Free replacements are warranted for a period of SIX(6)MONTHS.
- 3) Warranty does not cover the follow damage(s): product been modified, improper installation or use, improper clean electrical surges or nature caused damage.
- 4) Warranty does not cover the loss of third party with direct or indirect damage(s).

**Lifelab.**  
Longer Better Living

**Company: Lifelab BV**

 E-mail: [info@lifelab.shop](mailto:info@lifelab.shop)

 Address: Kerstraat 107a, 5411CK Zeeland, Netherlands

**For Super600/Super1200**

# Instruction for Super600/Super1200


Dear respected customer,

Thanks for choosing Super600/Super1200. We sincerely hope you are satisfied with our products and services. Please enjoy your lighting therapy according to the instruction below.

Please wear goggles or close eyes when enjoying therapy. Thank you very much.




## Operation gestures

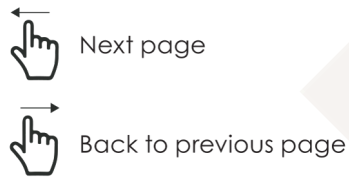

## Buttons



## Menu introduction

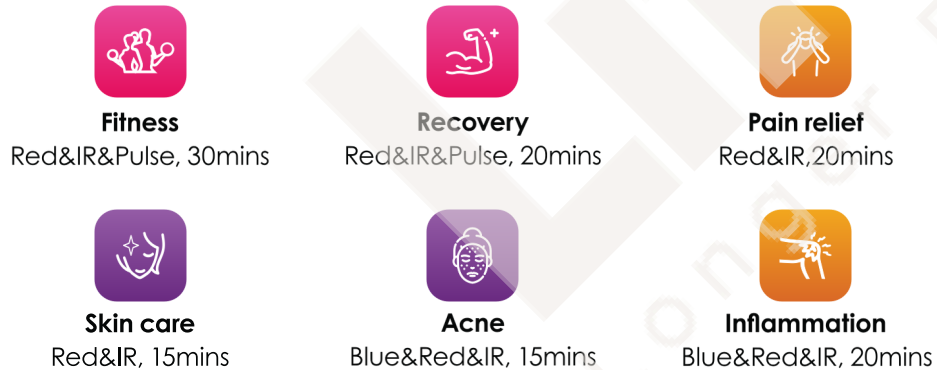


## Steps

**Step 1:**  
Unbox your device and accessories.



**Step 2:**  
Install the floor stand onto the panel.

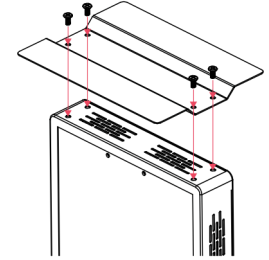


## Getting started

Plug in the device and switch it on, the screen will turn on.

**Default modes:**  
Choose your mode, then start the corresponding therapy.



**Customized mode:**

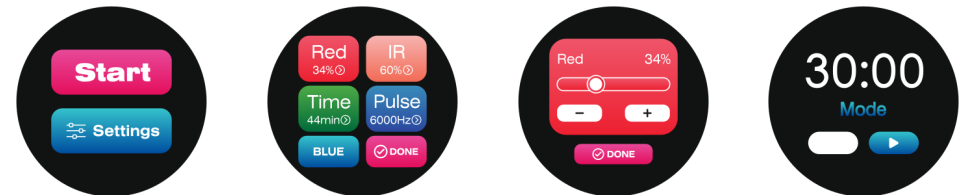

- Press "start" to begin the saved settings, first time with default settings (R100% / IR100% / no Blue / 20min / 0Hz)

- Press "settings" to set your values and press "done" to save. "Blue" is for with or without blue light.