





Instruction of Remote Controller

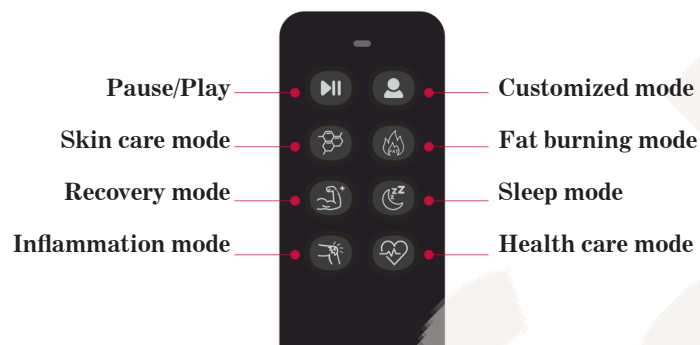

Before using the remote controller, please check if battery is included.

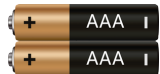

Battery normally not included for shipping reason

Warranty Clauses

- 1) Professional series are warranted against defects in materials and workmanship for a period of ONE(1) YEARS from the purchase date with free replacement.
- 2) Free replacements are warranted for a period of SIX(6)MONTHS.
- 3) Warranty does not cover the follow damage(s): product been modified, improper installation or use, improper clean electrical surges or nature caused damage.

Lifelab.
Longer Better Living

Company: Lifelab BV

✉ E-mail: sales@lifelab.shop

📍 Address: Kerstraat 107a, 5411CK Zeeland, Netherlands


For Professional Series

Instruction for Professional Series


Dear respected customer,
Thanks for choosing Lifelab BV. We sincerely hope you are satisfied with our products and services. Please enjoy your lighting therapy according to the instruction below. Please wear goggles or close eyes when enjoying therapy. Thank you very much.



Professional360

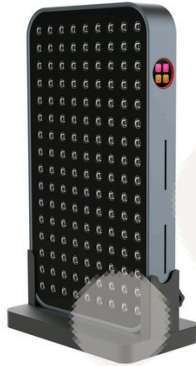

Professional720

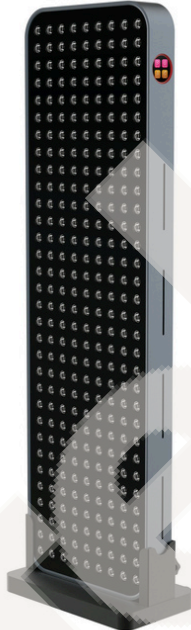

Professional1440

Operation Gestures



Next page



Back to previous page

Buttons



Play



Pause



Cancel

Menu Introduction



Skin care

Red&NIR&Green, 20mins
(10 mins R/NIR/G + 10 mins R/NIR)



Fat burning

Red&NIR&Pulse, 30mins



Recovery

Red&NIR&Pulse, 30mins



Sleep

Red&NIR, 20mins



Inflammation

Inflammation
Red&NIR&Blue&Pulse, 30mins
(15 mins R/NIR/B + 15 mins R/NIR)



Health care

Red&NIR, 30mins



Customized mode

Set customized values of Red/NIR/Blue/Green/Time/Pulse and save for quick use.

Getting Started

Plug in the device and switch it on, the screen will turn on.

Default modes:

Choose your mode, then start the corresponding therapy.



Customized mode:

- Press “start” to begin the saved settings.

- Press “settings” to set your values of spectra, pulse, time, and press “done” to save.

