




Voice Control

User Manual

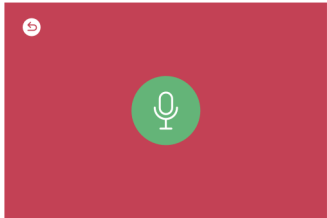

Turn on/off the voice control function

Voice settings

Wake-up: Hi/Hello iBeauty

- C: Run Skin care mode
- C: Run Acne mode
- C: Run Fat burning mode
- C: Run Hair Care mode
- C: Run Relax mode
- C: Run Sleep mode
- C: Run mode 1/2/3
- C: Stop
- C: Volume up/down
- C: turn on/off Bluetooth

- R: **Yes, I'm here. (Valid for 15s)**
- R: OK, Skin care mode is on
- R: OK, Acne mode is on
- R: OK, Fat burning mode is on
- R: OK, Hair Care mode is on
- R: OK, Relax mode is on
- R: OK, Sleep mode is on
- R: OK, mode 1/2/3 is on
- R: OK
- R: OK (10% each order)
- R: OK

Notes:

- 1) Make sure to consult with a physician before using the device if you have any of these conditions: acute diseases, febrile illnesses, blood pressure disorders, infectious diseases, history of skin cancer, skin disorders caused by diabetes or other diseases, currently pregnant or breast-feeding, sensitive skin due to allergic reactions, etc.,.
- 2) Please do not look directly at the device as it may cause discomfort or irritation to the eyes.

Warranty clause

- 1) Lifelab BY Ultimate medical series are warranted against defects in materials and workmanship for a period of ONE (1) YEAR from the purchase date with free replacement.
- 2) Free replacements are warranted for a period of SIX(6) MONTHS.
- 3) Warranty does not cover the follow damage(s): product been modified, improper installation or use, improper clean, electrical surges, or nature caused damage.
- 4) Warranty does not cover the loss of third party with direct or indirect damage(s)



Company: Lifelab BV

E-mail: info@lifelab.shop

Address: Kerstraat 107a, 5411CK Zeeland, Netherlands

For RL-Ultimate

Dear respected customer,
Thanks for choosing Lifelab BV. We sincerely hope you are satisfied with our products and services. Please enjoy your lighting therapy according to the instruction below. Please wear goggles or close eyes when enjoying therapy. Thank you very much.

Instruction for RL-Ultimate



Menu introduction

Skin Care
Red&IR, 30mins

Acne
Blue&Red&IR, 15mins

Fat Burning
Red&IR, 30mins

Relax
Red&IR, 20mins

Hair Care
Red&IR, 20mins

Sleep
Red&IR, 20mins

Customized mode 1, 2, 3
Set customized values of Red/IR/Blue/Time/Pulse and save for quick use.

Operation gestures

Next page

Back to previous page
Slide out the countdown interface from the main interface

Buttons

Back to menu

Play

Pause

Cancel

Getting started

Plug in the device and switch it on, the screen will turn on.

Default modes:

Choose your mode, set upper body / whole body / lower body, and begin therapy directly.

After the countdown interface returns, it can be slid out from the main interface using the gesture

Customized modes:

- Choose the mode, set upper body / whole body / lower body.
- Press " start" to begin the saved settings, first time with default settings (R100% / IR100% / no Blue / 20min / 0Hz)
- Press "settings" to set your values and press"done" to save. "Blue" is for with or without blue light.

Ps: you can preset all customized modes 1, 2, 3 for quick use.

Settings

Language

Brightness

Bluetooth

Turn on Bluetooth, find device"Beauty" in your phone, connect to play music.

Volume

Slide to adjust the volume